

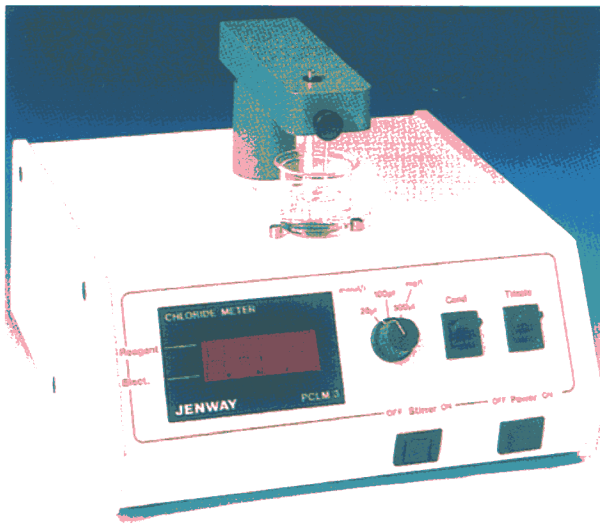
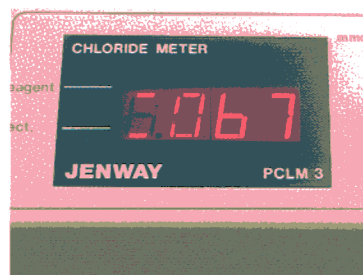


MODEL PCLM3 CHLORIDE METER

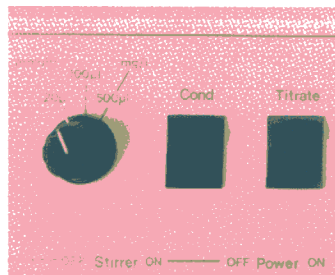

- Clinical & Industrial versions
- Three measurement ranges
- Large, clear LED display
- Simple to operate
- Fast, reliable results
- Low maintenance

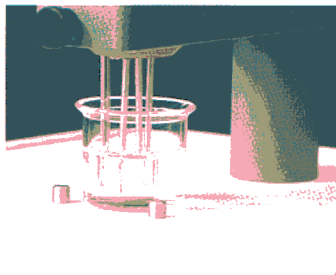

The Jenway PCLM3 is designed as a dual purpose chloride analyser for clinical and industrial applications. The unit has three measurement ranges; two are for use in clinical determinations and require $20\mu\text{l}$ or $100\mu\text{l}$ samples, and the third range is for industrial determinations and requires a $500\mu\text{l}$ sample.

The system operates on the proven coulometric principle, providing a quick and easy method to multiple measurements. Automatic indication is given when the reagent requires replacing or the electrodes need cleaning.

The PCLM3 will determine concentrations in biological samples such as serum, plasma, urine, sweat and CSF, and in industrial samples with low ionic strength, neutral pH and free of silver halides, other silver reactive components (other than chloride), low solid matter and low levels of dissolved solids.

ALYS Labware, Lausanne

Tel : 021 312 42 60

Fax : 021 312 42 61

labware@alys-technologies.com

Methodologies are available for the measurement of chloride in:

ALYS Labware, Lausanne

1. Oil based mud
2. Sodium bicarbonate crystals
3. Smoked fish
4. Biological fluids
5. Butter
6. Milk
7. Cement
8. Cheese
9. Curry powder & spices
10. Seasonings


Tel : 021 312 42 60

Fax : 021 312 42 61

labware@alys-technologies.com

ACCESSORIES

024 018	Glass beaker
025 011	Acid buffer solution (500ml)
025 012	Gelatin (30ml)
025 013	Chloride standard (clinical) (180ml)
025 014	Chloride standard (industrial) (180ml)
060 028	Electrode polish
060 029	Stirrer bar (qty 10)
504 052	Silver electrodes (1 cathode/2 detectors)
504 053	Silver anodes (qty 3)

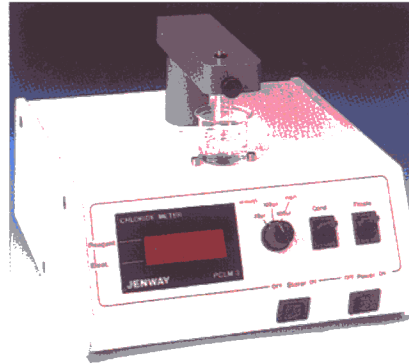

SPECIFICATIONS

	Clinical	Industrial
Range	10 to 299mmol/l	10 to 999mg/l
Sample volume	20 μ l or 100 μ l	500 μ l
Reproducibility	\pm 1% for 100 μ l sample @ 100mmol/l \pm 1.5% for 20 μ l sample @ 100mmol/l	\pm 3mg/l
Linearity	Better than \pm 1mmol/l or \pm 1% of the concentration value (whichever is the greater) over the range of 10 to 299mmol/l	\pm 3mg/l or 1% (whichever is the greater) over the range of 10 to 999mg/l
Power	240/115V a.c. \pm 10% @ 50/60Hz	240/115V a.c. \pm 10% @ 50/60Hz
Size	240 x 215 x 160mm	240 x 215 x 160mm
Weight	3.1Kgs	3.1Kgs
504 043	Model PCLM3 (clinical version) supplied complete with mains lead, buffer, gelatin, clinical chloride standard, 3 silver electrodes, 2 silver anodes, electrode polish, glass beaker, stirrers, tweezers and manual.	
504 044	Model PCLM3 (industrial version) supplied complete with mains lead, buffer, gelatin, industrial chloride standard, 3 silver electrodes, 2 silver anodes, electrode polish, glass beaker, stirrers, tweezers and manual.	

UK & Europe
Jenway Limited
Gransmore Green
Felsted, Dunmow
Essex CM6 3LB
ENGLAND
tel: 01371 820122
fax: 01371 821083
e-mail: sales@jenway.com

North & South America
Techne Incorporated
743 Alexander Road
Princeton, NJ 08540-632
UNITED STATES
tel: (609) 452-9275
fax: (609) 987-8177
toll free: 800-225-9243

Jenway Middle East Limited
P.O. Box 27842,
Nicosia
CYPRUS
tel: 357 2 660423
fax: 357 2 660424